

# MicroBil

## Bilirubin Analyser



- ? BILIRUBIN Analyser MicroBil of Krish Biomedical provides simple, Quick and accurate results of Total Bilirubin.
- ? MicroBil is advanced and optimised instrument for measuring Neonate's Bilirubin.
- ? Microbil eliminates interference from turbidity and hemolysis with its unique dual wavelength measurement system.

**Without any reagent  One Drop of blood is required for test ,**


- ?Microbil requires minimal sample volume
- ?Built in self Diagnostics
- ?Life long LED light source that does not require replacement
- ?New Design capillary holder that prevents damage & breakage of Capillary tubes.
- ?Quick Reading photometer displays results in 2 seconds.
- ?Built in Battery Backup that can last for over 500 Test.

### SPECIFICATIONS:

Light Source	Life Long Light Emitting Diodes
Wave length	460nm & 570 nm
Measuring Range	0-40 mg/dL
Display	Alpha numeric Liquid Crystal Display
Accuracy	Within $\pm 3\%$
Power Supply	12V 1A Remote Power Supply
Dimension & Weight	222(W) x 180(D) x 100(H) mm & 1.5 Kg
Built in Battery Backup	Sufficient for 500 Test

### STANDARD ACCESSORIES:

Remote Power Supply	1 pc
Capillary Tube, Heparinized	1 Box (400 pcs)
Sealing Wax	1 pc
Capillary Adaptor	3 pcs
User Guide	1 pc

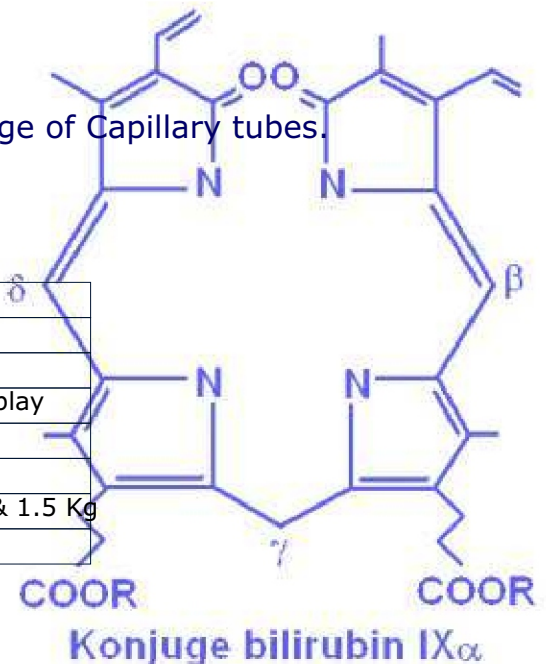

# Krish Biomedicals

**An ISO 9001:2008 Company**

**Address : T-7, Pankaj Central Plaza, Near Plato Public School,  
Patparganj, New Delhi - 110092**